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

EN
USER INSTRUCTIONS
Refer to the product label/markings for detailed information on the corresponding standards.
Only standards and icons that appear on both the product and the user information below are applicable. All these products comply with the requirements of Regulation (EU) 2016/425 and Regulation 2016/425 as brought into UK law and amended and compliant with standard EN149:2001+A1:2009.

1. This filtering device is tested according to EN149:2001 and performs the requirements of its corresponding filter- class. The mask protects against solid and oil- based aerosols, but not against solvents for spray painting operations, carcinogenic and radionuclide substances, micro-organism, viruses and substances with biochemical potential, gases and vapour with harmful effects.
2. No alteration or modification is allowed. Never use damaged masks, in that case discard these masks.
3. The wearer of the mask must be well trained to use properly before starting work.
Warning:
Facial hair (beard or sideburns) of the user within the area of the sealing edges of the mask may jeopardize the protection efficiency of the mask itself.
4. The safe functioning of this filtering device depends on the ambient atmosphere. This atmosphere must contain at least 19.5% oxygen per volume, therefore it is not allowed to use this filtering device in poorly ventilated areas or wherever a deficiency of oxygen may occur. Do not use this mask, when the concentration of contaminants is immediately dangerous to life or health. The concentration, kind of substance(s) and harmfulness characteristics of the substances in atmosphere must be known. Leave the dangerous area immediately, if you feel distressed or experience difficulty in breathing.
5. **Warning:**
The device is marked "R D" meaning it's a reusable filtering half mask**.** It should be discarded and replaced when the breathing resistance becomes excessive, the cleaning and disinfecting is not complete, or if it becomes deformed after used. The masks should be stored in a clean, airtight container after each use or after cleaning or disinfecting.
6. These instructions for use are a guide to the correct use of the product and to prevent health risks. They must be read, completely understood and carefully followed, in the same way as all national laws or other regulations in forces regarding the prevention of accidents at work. The filtering devices must be stored and transported in their original package and protected by the storage temperature and humidity as suggested by the manufacturer.
CLEANING AND DISINFECTING
The device is a reusable mask, so it can be used for more than one shift.It must be cleaned at the end of a shift by following steps:
1. Use one hand to hold the mask horizontally above a desk with the outside facing down.
2. Tap the mask against the desktop 10 times.
3. Use 75% commercial alcohol pad (no-drip humidity) to wipe clean the valve surface and wait until dry before reusing.
** Do not immerse product in water
** Store cleaned respirator in clean, airtight container.
WARNINGS
1. Warning: Do not use this mask as protection against harmful gases or vapours, solvents from spray-painting operations, or in explosive atmospheres, or in an environment where harmful gas is present.
2. Do not use this mask in tanks, enclosed places or in any circumstances where harmful gas is likely to be present or in oxygen deficient atmospheres (< 19.5% in volume.)
3. This mask should not be worn when conditions prevent a good face to facepiece seal. Such conditions include the growth of beards or moustaches and sideburns and the use of spectacles or any other devices containing temples or straps which pass between the mask sealing edge and the wearers face.
4. Leave the work area immediately if breathing becomes difficult or if dizziness or any other signs of distress occur.
5. The instructions for use and warning are helpful for correctly using the product and to prevent risks to health. It is peremptory to read and follow them carefully.
6. Checks must be done prior to use the mask. i.e.: dirty, damaged, expiry date etc.
FITTING INSTRUCTIONS (SEE DIAGRAM - D1)
Before fitting the respirator, carefully inspect the entire respirator to ensure that the respirator is not damaged and is functioning correctly. Discard the entire respirator if any part does not function properly.
1.Place the respirator over your nose and mouth. Be sure the nose shape position is on top.
2.Pull the bottom strap over your head until it rests below the ears at the back of your neck. Untwist the strap.
3.Pull the top strap over your head until it rests on the crown of your head above your ears. Untwist the strap.
4.Without removing the respirator, adjust the straps until the respirator is secure. To tighten the respirator, gently pull the ends of both straps. To loosen the respirator, press on the back side of the buckle.
5. To check the correct fit, cup the respirator with both hands and inhale sharply.
If any airflow indicates leakage around the edges of the respirator, stretch the headbands again by pulling them backwards. Then correct the position of the respirator (see picture) and check again. WARNING! It is not allowed to use the mask without proper fit.DO NOT ENTER A CONTAMINATED AREA WITHOUT A PROPER SEAL.
** Note: Checking the nose shape position using one hand may result in improper fit and less effective respirator performance. Use two hands.
STORAGE AND SHELF LIFE
The shelf life of the mask is 3 years and the mask should be placed in a container and stored in a clean cool and dry atmosphere to ensure that it is protected again dust, chemicals, moisture, excessive heat and direct sunlight.

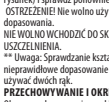TEMPERATURE RANGE OF STORAGE CONDITIONS

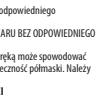MAXIMUM RELATIVE HUMIDITY OF STORAGE CONDITIONS

Download declaration of conformity @ www.portwest.com/declarations

MANUFACTURER
PORTWEST

Portwest UC, IDA Business Park Westport, Co Mayo, F28 FY88, Ireland

NAME & ADDRESS OF THE NOTIFIED BODY HAVING ISSUED EU CERTIFICATE:
INSPEC INTERNATIONAL B.V. - NOTIFIED BODY: 2849
BEECHAVENUE 54-62, 1119 PW, SCHIPHOL-RIJK, NETHERLANDS
Model : P306

NAME & ADDRESS OF THE APPROVED BODY HAVING ISSUED UKCA CERTIFICATE:
INSPEC INTERNATIONAL LTD, NO. 0194
56 LESLIE HOUGH WAY, SALFORD, GREATER MANCHESTER, M6 6AJ, UK
Model : P306

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CE 2849 UKCA 0194

EN149:2001+A1:2009

See Information supplied by manufacturer

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